



USB Heating Eye Mask **iRelief™ VEM-200**

Provide Effective Warm Compress to Dry Eyes

Listen to the clinical

Product Concept

iRelief™ is a smart eye mask designed to provide heating comfort to the eyes. It works by connecting eye mask with power bank or computer or any other power source. It can relieve eye fatigue and discomfort caused by dry eye or other eye problems. It can also accelerate the blood circulation around the eye to improve black eye circle conditions. There is a removable sheet on the eye mask for users to clean and keep hygiene. The sheet is machine washable and dryable.



Temperature
Adjustable



Nasal Arch Fit



Machine Washable




Auto Switch off



FlexStrip™ Fits Face Better

FlexStrip™ is innovatively designed and applied in the nose bridge area of the iRelief™ eye mask. FlexStrip™ is adjustable according to the nose shape to make the iRelief™ eye mask better fit on the face. Better fit provides better heating effect at the eye lid area.

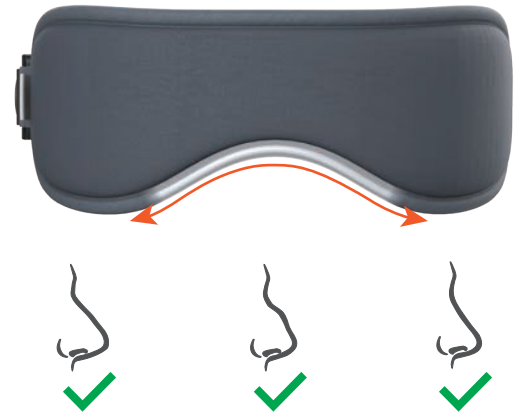


- 39°C
- 41°C
- 44°C

3 Heating Options with Smart Control

iRelief™ provides 3 optional heating temperatures to meet different comfort demands. The heating time is 10 minutes for each option. When the 10 minutes' heating SPA is finished, it will automatically switch off.

High Quality Material in Each Layer

Each layer of the eye mask is made up with high quality materials to ensure comfortable and safe wearing experience.

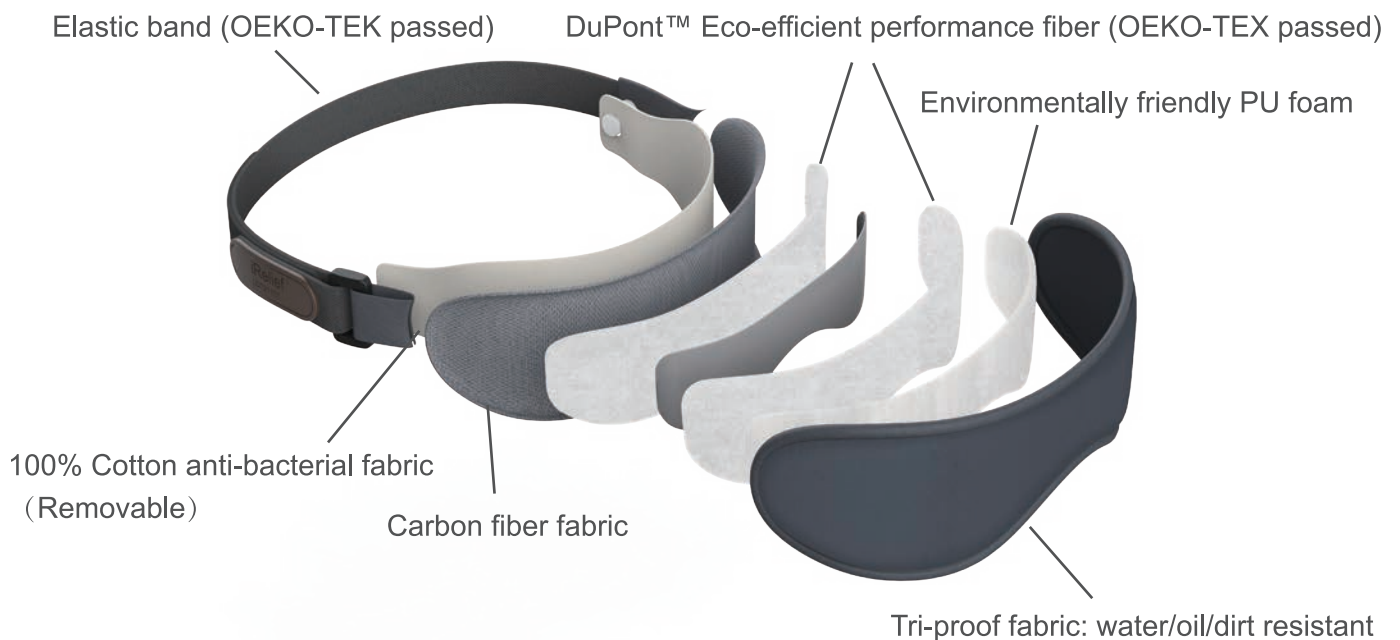

Specification

VEM-200

Heating method	Powered by USB connection
Heating temperature	39°C, 41°C, 44°C
Heating duration	12 minutes
Input voltage	DC 5V/1A
Power	4W (max.)
Color	Space grey
size	22cm x 7cm
Weight	25g

DuPont™ is the trademark of DuPont de Nemours, Inc