

NEO COSMO COLOR SOFT CONTACT LENS



NEO VISION
COLOR SOFT CONTACT LENS



NEO VISION

COLOR SOFT CONTACT LENS



CEO Greetings;

Our Innovative products developed based on our long experience in the research of the new various contact lenses technologies have been introduced to the markets, and we are always making our repeated hard work to supply our customers with the good quality products. They have already achieved an excellent reputation for their high quality in the European and Japanese markets. We will do our best to grow up to be one of the most reliable leading firms that is manufacturing contact lenses in the world, and run a well based company with the positive management and good quality control system to satisfy our customers needs.

Thank you.

Sincerely yours
Neo Vision Co.,Ltd.

Kyung Hwa Kim / President



CIRTIFICATE



CE certificate



ISO certificate



SDA certificate



certificate of KFDA



certificate of INNO-BIZ

- 1994. 07. To establish Poong Gwang Optical Research Laboratory by the assistance of high polymers lesson, technology of Jeon Nam University.
- 1999. 11. Since beginning to export to Singapore, we could achieve a great work of the diversification of export markets.
- 2001. 09. To move the office and change the name of the company to Neo Vision.
- 2002. 10. To be awarded the certificate of SDA for sales in China and to establish an individual cooperation.
- 2003. 07 EN ISO 9001 certified. (UNK20559)
- 2004. 07. To secure CE certificate. (EN ISO 13488:2000, EN46002:1996 etc.)
- 2005. 01. To share technology of medical ozonizer with Grow Co. in Japan.
- 2005. 04. To organize the committee of academic-industrial cooperation with department of optical science in Yeo Joo University.
- 2005. 05. To apply for a patent of cosmetic color lens material (License No.10-2005-004792) and coating technology.
- 2005. 08. To regist as venture company (Business Administration) Reg. No. 51627032-1-00992.
- 2005. 09. To establish a R&D center.(NO. 20052784)
- 2006. 04. To apply again for articles of incoporation and business registration. (To add medical device).
- 2006. 06. To start processing on commission with Chinese factory.
- 2006. 09. To get supply of academic-industrial cooperation with department of optical science in Chodang University.
- 2006. 10. To apply for a patent of design of lenses. (License No. 30-0249476)
- 2006. 11. To be awarded INNO-BIZ (No. 6061-2883) by Business Administration.
- 2006. 11. To obtain a patent for Sandwich process. (License No. 10-0647133)
- 2007. 02. To be awarded KGMP. (No. MGK-25)
- 2007. 11. To be awarded Small and Medium enterprises by Small and Medium Business Administration.
- 2008. 04. To establish the second factory with cast mold system in Wonju city.


NEO COSMO MONTHLY

DISPOSABLE COLOR CONTACT LENS



BLISTER PACK-6PCS

3 CURVES

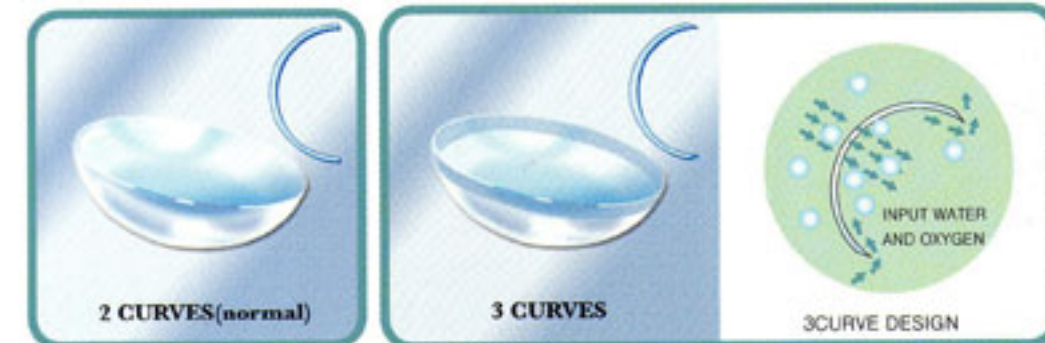

An epoch-making aspheric design makes tearing more in a smooth cycle and the external design has been made out with 3curves to help the patients have the most comfortable feeling.

S.P Process



S.P Process is to infiltrate the cosmetic color between the inside and outside layer of the contact lens, which makes the wearing feeling comfortable and the color is not decolorized.

NEO COSMO MONTHLY SPECIFICATION

B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	45%
Color	Blue, Green, Violet, Brown, Gray Aqua, Hazel, Honey
Manufacturing techniques	casting mold

EYEPIA MONTHLY


THE FEATURES OF CONTACT LENSES
RECENTLY UPGRADED




EYEPIA MONTHLY SPECIFICATION	
B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	55%
Color	BlueTinted
Manufacturing techniques	casting mold

TORIC LENS

FRONT TORIC DYNAMIC

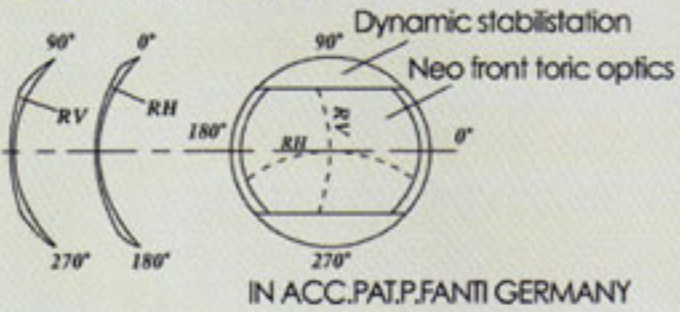

BACK TORIC PRISM

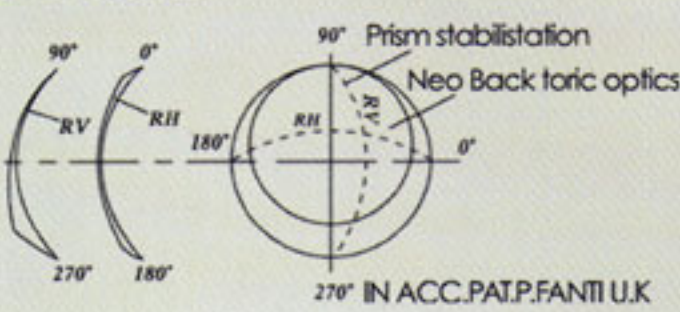

FRONT TORIC LENS has no problem in the movement of axis and it has the feature which enables to promote the stability of axis during wearing the lenses. Dynamic front toric lens is auto free focus lens and corneal trouble has been removed when using full automatic process by the computerized design.

BACK TORIC LENS has been manufactured with semi casting mold system and it's axis of weight center is placed in the lower position. Toric lens gives comfortable feeling for beginners to wear, because of no movement of axis.

FRONT, BACK TORIC SPECIFICATION	
B.C(mm)	8.6
DIA	14.0
Power	0.00 ~ ±10.00
Cycle	-0.75~0.35
Axis	0 ~ 90 ~ 180
Materials	2-HEMA
Water contents	38%,42%
Color	Blue Tinted
Manufacturing techniques	Semi casting mold

SPHERICAL LENS


FOUR-TONE COLOR

THE COSMETIC LENS THAT NO ONE CAN IMITATE



NEO COSMO FOUR-TONE
COLOR LENS It has taken two
years to complete it since we
started to develop this 4T lens
in our laboratory.
The design of lens is entirely
composed of the light and
dark tones for each color
so that it can naturally show all
the details of shapes and colors.

NEO COSMO COLOR SOFT CONTACT LENS

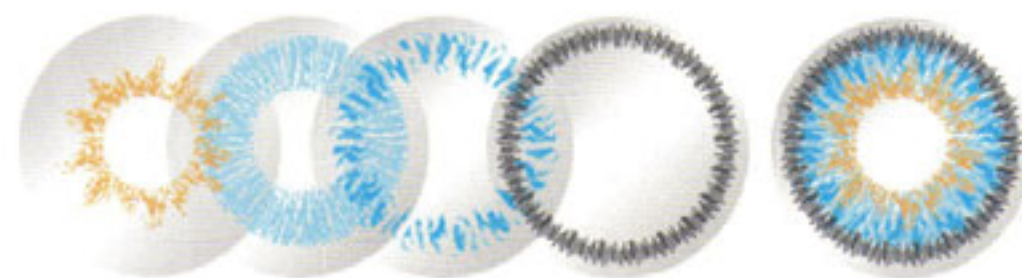

4 COLOR SYSTEM



BLUE



N411



GREEN



N412



BROWN



N414



GRAY



N415

NEW



BLUE



N421



GREEN



N422



BROWN



N424



GRAY



N425



4 COLOR SYSTEM

FOUR-TONE COLOR SPECIFICATION

B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	45%
Color	Blue, Green, Brown, Gray
Manufacturing techniques	Semi Casting mold,Casting mold



FOUR-TONE OF COLOR CONTACT LENS

TRI-TONE COLOR

REPRESENT YOUR IDEA

WITH THE CHARMING EYES



WWW.NEOLENS.CO.KR

NEO COSMO TRI-TONE COLOR is a natural color lens based on the light brown pattern of 3 tone inside diameter, Which is suitable for other colors. The color of eyes during wearing the lens also shows the natural beauty to others, and it is a new concept lens of NEO VISION which doesn't give the eyestrain at all, providing the patient with the comfortable feeling to wear.

NEO COSMO COLOR SOFT CONTACT LENS



3 COLOR SYSTEM



BLUE



N331



GREEN



N332



VIOLET



N333



BROWN



N334



GRAY



N335



AQUA



N336



HAZEL



N337



HONEY



N338



BLUE



N341



GREEN



N342



VIOLET



N343



BROWN



N344



GRAY



N345

TRI-TONE COLOR SPECIFICATION

B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	45%
Color	Blue, Green, Violet, Brown, Gray, Aqua, Hazel, Honey
Manufacturing techniques	Semi Casting mold, Casting mold



TRI-TONE COLOR

TWO-TONE COLOR

THE BEAUTIFUL EYES WITH
CIRCLE AND FLOWER PATTERN



WWW.NEOLENS.CO.KR

NEO COSMO TWO-TONE
COLOR LENS is well-matched
with circle and flower pattern.
It must be the most suitable
lens for the young people
who want new design and
clear eye look.
Cosmetic lens with circle
concept make your eyes
looking larger than actual.

NEO COSMO COLOR SOFT CONTACT LENS



2 COLOR SYSTEM



N241

BLUE



N244

BROWN



N245

GRAY



N221

BLUE



N222

GREEN



N223

VIOLET



N224

BROWN



N225

GRAY



N226

AQUA



N227

HAZEL



N228

HONEY

TWO-TONE COLOR SPECIFICATION

B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	45%
Color	Blue, Green, Violet, Brown, Gray, Aqua, Hazel, Honey
Manufacturing techniques	Semi Casting mold, Casting mold



TWO-TONE COLOR

CIRCLE & EMOTION

CIRCLE LENS AND EMOTION LENS
FOR NEW GENERATION



WWW.NEOLENS.CO.KR

CIRCLE LENS make your eyes bigger and clearer.
EMOTION LENS is specially designed with symbolizing the shiny tears with four different colors. Natural colors show your beautiful emotion to others.



RING
CIRCLE



N010



BLACK
CIRCLE



N011



CHAGALL
CIRCLE



N012



DALI
CIRCLE



N013



EMOTION
GRAY



N026



EMOTION
BROWN



N027



EMOTION
BLUE



N028



EMOTION
PINK



N029



N021



N022



N023



N024



N025



CARRA
CIRCLE



N014

CIRCLE & EMOTION SPECIFICATION

B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	45%
Color	Blue, Green, Violet, Brown, Gray Black, Brown Black, Pink
Manufacturing techniques	Semi Casting mold, Casting mold



THESE LENSES ARE SPECIALLY DESIGNED WITH ULTRA SINK MOLDING PROCESS

TWO DOT COLOR



N231



N232



N233



N234



N235

TRI COLOR



N311



N312



N313



N314



N315

TRI DOT COLOR



N321



N322



N323



N324



N325

NEO COSMO SPECIFICATION

B.C(mm)	8.6
DIA	14.2
Power	0.00 ~ ±10.00
Materials	2-HEMA
Water contents	45%
Color	Blue, Green, Violet, Brown, Gray
Manufacturing techniques	Semi Casting mold, Casting mold



WWW.NEOLENS.CO.KR

NEO PLUS MULTI & DROP

NEO PLUS MULTI

(145ml/360ml)

NEO-PLUS, multi-purpose contact lens solution has the functions such as; disinfection, cleaning & rinsing, protein removal and preservation which don't have any sore eyes when in use.

EYE SENSE DROP

Eye-sense drop is an effective lubricant that provides the long-lasting wearing for wearer feeling comfortable in the eyes.

Product development & Production

A. Product range


1. Contact lens - Moisturing Lens & Protein Lens
2. General Lens Care
 - Saline
 - Multi Solution
 - Lens Drop (lubricant eye drop)

B. Study & Research

1. Medical Supplies (eyewash)
 - Drying Eye
 - Cataract / Congenital [acquired] Cataract
 - Glaucoma
2. Bio Field - Drug Delivery System
 - Artificial Iris
 - Characteristic Contact lens



• NEO PLUS 4 STEPS SYSTEM

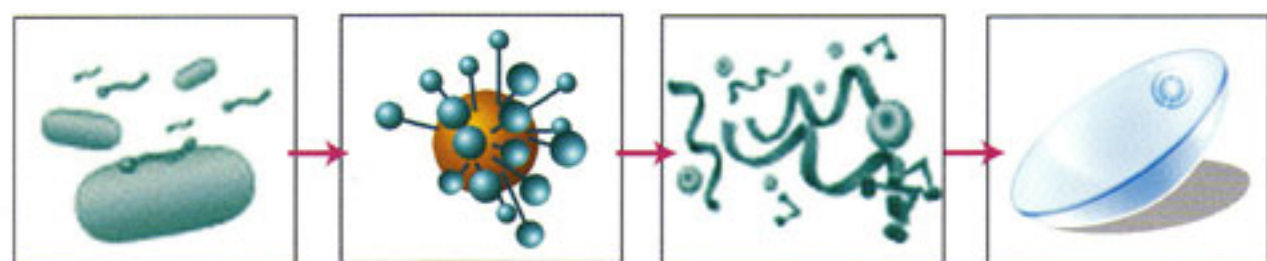

- Step 1 : Disinfecting polyhexamide
- Step 2 : Cleaning Poloxamer 407
- Step 3 : Removing proteins Poloxamer 407 EDTA
- Step 4 : Preserving-potassium chloride, boric acid buffer

• ANTIPOLLUTION CAP



- Step 3 : Clean and rinse Drop NEO PLUS solution on the lens, and clean it up
- Step 1 : safety cap Turn the safety cap to the direction to open.
- Step 2 : Antipollution cap Turn the antipollution cap to the direction to open

SUTRA
OPEN

NEO VISION



P.O. Box 13893
Laudium 0037
Pretoria, South Africa

Tel: (012) 370 4175/ (012) 370 3383
Fax: (012) 370 1557
Cell: 082 414 1472
sales@eurotechoptical.com
www.eurotechoptical.com

Exclusive distributors of
Optical / Ophthalmology equipment-
Autorefractor / Keratometer, Slitlamps
Fundus Cameras and other optical/
ophthalmology products at competitive
prices. Contact lenses, colour toric
lenses & optical frames available.

Contact us for the most competitive prices