


Distance Acuity Wall Charts

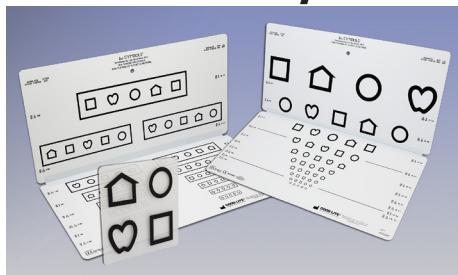


LEA SYMBOLS® 10-Line Folding Distance Chart (10 feet/3 meters)

This eye chart folds in half for added mobility or easy storage.

Proportionally-spaced lines; line sizes range from 20/200 to 20/16 (6/60 to 6/4.8) equivalent. The 20/160 (6/48) and 20/100 (6/30) lines have been omitted. 50% spaced rectangle (Massachusetts Visual Acuity Test format) with line sizes from 20/80 to 20/16 (6/24 to 6/4.8) equivalent. Has a hole for hanging the chart on the wall. Includes response key and instructions. 16" x 16" (40.6 cm x 40.6 cm). Folds to 16" x 8" (40.6 cm x 20.3 cm).


..... 250200

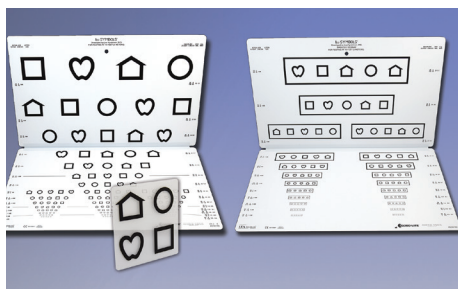

LEA SYMBOLS® 15-Line Distance Chart (10 feet/3 meters)

This test is designed for vision screening at schools and assessment of vision of children with special needs beginning with age of 3 years. Proportionally spaced lines; line sizes range from 20/200 to 20/8 (6/60 to 6/2.4) equivalent. 50% spaced rectangle (Massachusetts Visual Acuity Test format) with line sizes from 20/100 to 20/8 (6/30 to 6/2.4) equivalent. Has a hole for hanging the chart on the wall. Includes response key and instructions. 17" x 21" (43 cm x 53.4 cm)..


..... 250100

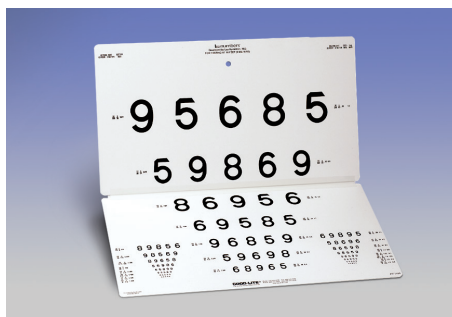

LEA NUMBERS® 15-Line Distance Chart (10 feet/3 meters)

Proportionally spaced (logMAR) lines; line sizes range from 20/200 to 20/8 (6/60 to 6/2.4) equivalent, 0.10 to 2.50. Has hole for hanging on the wall. 18" x 20" (45.7 cm x 50.8 cm). Folds to 18" x 10" (45.7 cm x 25.4 cm).


..... 271100

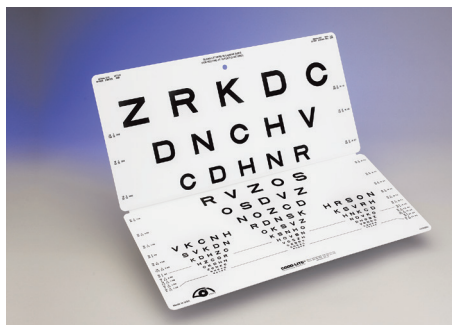

Sloan Letters Folding Chart (10 feet, 3 meters)

Proportionally spaced (logMAR) lines; line sizes range from 20/200 to 20/8 (6/60 to 6/2.4) equivalent, 0.10 to 2.50. 18" x 18" (45.7 cm x 45.7 cm); folds to 18" x 9" (45.7 cm x 22.9 cm).

..... 735000

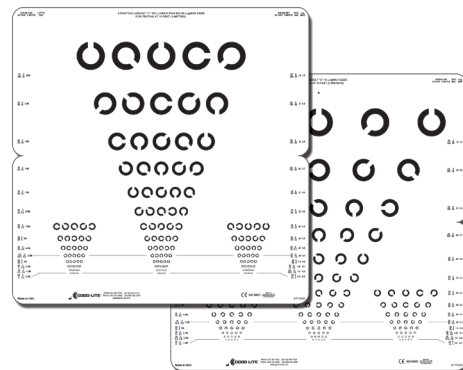

8 Position Landolt C Folding Wall Chart (10 feet/3 meters)

8 Position Landolt C chart with C openings at up, down, left, right and each 45° angle. 8 Position Landolt C wall chart is printed on both sides. The front contains proportionally spaced (geometric progression) lines while the backside has 25% linear spaced optotypes. The chart is perfect for vision testing those who maybe illiterate or unfamiliar with letter or number optotypes.

Both sides have line sizes ranging from 20/100 to 20/8 (6/30 to 6/2.4) equivalent. Hole at top of chart for hanging on the wall.

Chart measures 18" x 18" (45.7 cm x 45.7 cm) and folds in half to 18" x 9" (45.7 cm x 22.8 cm) for easy storage.

..... 779200



Cardiff Pediatric Acuity Test

Designed specifically for acuity measurement in children ages 6 months to 3 years. The test may also be helpful in other age groups for people with intellectual impairment, dementia, head injury, or stroke. Cards can be used at either 1 meter or 0.5 meter. This card set includes 36 cards providing 24 acuity levels when used at both distances and Children's Fun Frame Occluders.

Cardiff Pediatric Acuity Test 474000

Cardiff Low Vision Cards 475400


Cardiff Low Contrast Cards 566700

Insist On Genuine LEA Products by GOOD-LITE®

lea SYMBOLS®
lea-numbers®

Distance Visual Acuity Testing

Distance Acuity Handheld Tests

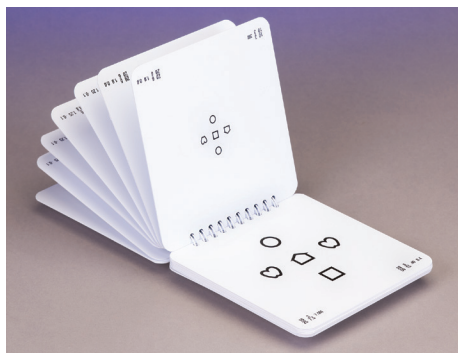

LEA SYMBOLS® Crowded Symbol Book (10 feet/3 meters)

Variation of the line test for children or adults who have difficulty concentrating on a chart. The center symbol is surrounded on four sides by symbols of the same size. The interaction between the center symbol and the surrounding symbols is approximately the same as in a line test. The test size is convenient and easy to administer during screenings, home visits, and unique testing situations. Often, this is the first line test that can measure visual acuity of a young child with amblyopia, because the test produces the effect of crowding. Line sizes range from 20/63 to 20/16 (6/19 to 6/4.8) equivalent, 0.32 to 1.25. Includes response key and instructions. 16 offset pages for easy flipping. 5" x 5" (12.7 cm x 12.7 cm).


..... 250700



LEA SYMBOLS® 50% Spaced Crowding Bars Book (10 feet/3 meters)

This book features 50% spaced contour interaction bars to provide single optotype acuity with the proper crowding effect. Optotype line sizes range from 20/100 to 20/20 (6/30 to 6/6) equivalent, 0.20 to 1.00. Includes response key and instructions. 17 offset pages for easy flipping. 4" x 4" (10.2 cm x 10.2 cm).


..... 259200

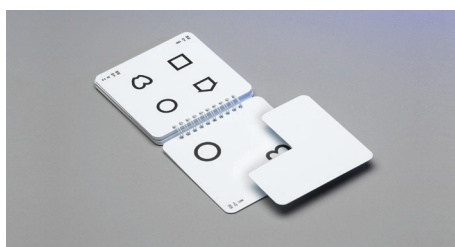

LEA SYMBOLS® Single Symbol Book (10 feet/3 meters)

The Single Symbol Book is ideal for people who cannot perform when tested with a line test. This book can help determine the best resolution of the visual pathway in a child/adult with amblyopia or impaired vision. Single optotype line sizes range from 20/200 to 20/8 (6/60 to 6/2.4) equivalent, 0.10 to 2.50. Includes isolating mask, response key and instructions. 13 offset pages for easy flipping. 5" x 5" (12.7 cm x 12.7 cm).

..... 250600

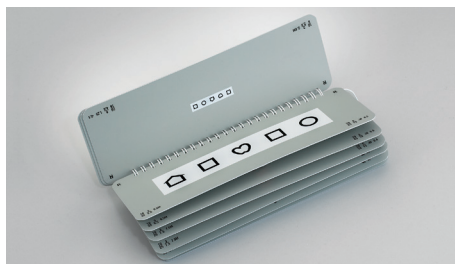

LEA SYMBOLS® Distance Vision Screener (10 feet/3 meters)

This test is designed for amblyopia testing or screening. The background is gray with the exception of a white rectangle spaced at 50% surrounding each line of symbols. Line sizes range from 20/80 to 20/16 (6/24 to 6/4.8) equivalent, 0.25 to 1.25. Includes response key and training cards. Pages are offset for easy flipping. Instructions are printed on the inside cover. 10 pages, 3.25" x 9" (8.25 cm x 22.9 cm).

..... 259900



LEA SYMBOLS® Flash Cards (10 feet/3 meters)

This set of flash cards was designed for testing distance visual acuity in children and adults who have difficulty with visual acuity charts, the Crowded Symbol Book test or the Single Symbol Book test. Line sizes range from 20/200 to 20/8 (6/60 to 6/2.4) equivalent, 0.10 to 2.50. Includes response key and instructions. 30 cards per set, 4" x 5" (10 cm x 12.7 cm) to 5" x 5" (12.7 cm x 12.7 cm).

..... 252700



Single HOTV Letter Book With 50% Spaced Crowding Bars (10 feet/3 meters)

This book features 50% contour interaction bars to provide single optotype acuity with the proper crowding effect. Optotype line sizes range from 20/100 to 20/16 (6/30 to 6/4.8) equivalent, 0.20 to 1.25. Includes response key, and instructions. 17 offset pages for easy flipping. 4" x 4" (10.2 cm x 10.2 cm).

..... 600350



LEA NUMBERS® Single Number Book (10 feet/3 meters)

The LEA NUMBERS® Single Number Book is ideal for situations in which the child/adult cannot perform when tested with a line test. Single optotype line sizes range from 20/200 to 20/8 (6/60 to 6/2.4) equivalent, 0.10 to 2.50. Includes isolating mask and instructions. 13 offset pages for easy flipping. 5" x 5" (12.7 cm x 12.7 cm).

..... 274000



Magnetic Sloan Letters Chart (20 feet/6 meters)

This durable chart is made of plastic-coated, galvanized steel with 2 sets of 20 magnetic strips (40 total) with 5 characters on each strip. Meets all qualifying standards for military aviators (e.g., officer candidates, pilots, air crewmen, etc). Approved by the Navy Bureau of Medicine and Surgery (BUMED). Chart measures 9.5" x 11" (24.1 cm x 28 cm).

20/20 size 400500
20/40 size (not shown) 400900